


Technical Data

REF 520

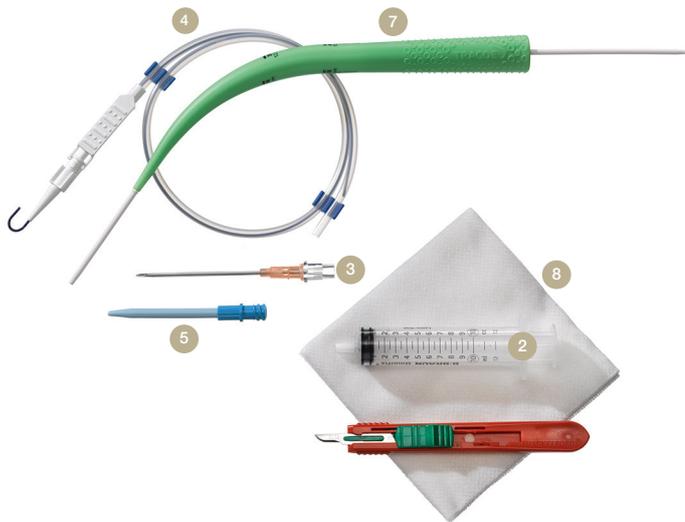

TRACOE experc Dilation Set for percutaneous tracheotomy

- 1 scalpel
- 2 syringe
- 3 puncture needle (14 G) with Teflon catheter
- 4 Seldinger guide wire made of kink-resistant nitinol with inserter
- 5 short dilator (14 Ch/Fr)
- 6 guiding catheter with safety stop
- 7 TRACOE experc dilator with hydrophilic coating (slides easily after moistening)
- 8 4 compresses

◆ in sterile package

The Innovation

The distinctive feature of the TRACOE experc set is the minimally traumatic inserter for the placement of the tracheostomy tube after percutaneous tracheostomy. There is a foldable silicone sleeve at the tip to smoothly bridge the gap in diameter between the inserter and the end of the tube.

Minimal traumatic tube placement

Inserting the tracheostomy tube via an inserter is a particular challenge of percutaneous dilation tracheotomy. Our primary tubes are supplied already pre-installed on the respective minimally traumatic inserter, which is adapted to the specific tube shape and length. After placing the tube in the trachea, the inserter is pulled out of the tube. The silicone sleeve flaps back so that the inserter can be removed.

The experc system offers the possibility to order the dilation set individually or in combination with a so-called P-tube from the product lines twist, twist plus, vario and vario XL as a set.

Technical Data

experc Sets *twist*



REF 320

TRACOE experc Set *twist* including REF 301-P with low-pressure cuff and minimally traumatic inserter



REF 321

TRACOE experc Set *twist* including REF 302-P with low-pressure cuff, fenestrated, and minimally traumatic inserter



REF 322

TRACOE experc Set *twist* extract including REF 306-P with low-pressure cuff, subglottic suction line and minimally traumatic inserter



REF 888-322

TRACOE experc Set *twist* extract including REF 888-306-P with low-pressure cuff, fenestration, subglottic suction line and minimally traumatic insertion system

To order, please specify: REF + Size, e.g. REF 320-08

experc Sets *twist plus*



REF 330

TRACOE experc Set *twist plus* including REF 311-P with low-pressure cuff and minimally traumatic insertion system



REF 331

TRACOE experc Set *twist plus* including REF 312-P with low-pressure cuff, double fenestrated, and minimally traumatic insertion system



REF 332

TRACOE experc Set *twist plus* extract including REF 316-P with low-pressure cuff, subglottic suction line and minimally traumatic insertion system



REF 888-332

TRACOE experc Set *twist plus* extract including REF 888-316-P with low-pressure cuff, double fenestrated, subglottic suction line and minimally traumatic insertion system

To order, please specify: REF + Size, e.g. REF 330-09

Technical Data

experc Sets *vario*



REF 420 

TRACOE experc Set vario including REF 450-P spiral-reinforced, with adjustable neck flange, low-pressure cuff and minimally traumatic insertion system



REF 422 

TRACOE experc Set vario extract including REF 470-P with adjustable neck flange, low-pressure cuff, subglottic suction line and minimally traumatic insertion system



REF 423 

TRACOE experc Set vario including REF 460-P with adjustable neck flange, low-pressure cuff and minimally traumatic insertion system

To order, please specify: **REF + Size**, e.g. **REF 420-08**

experc Sets *vario XL*



REF 421 

TRACOE experc Set vario XL including REF 451-P extra-long, spiral-reinforced, with adjustable neck flange, low-pressure cuff and minimally traumatic insertion system



REF 425 

TRACOE experc Set vario XL extract including REF 471-P extra-long, with adjustable neck flange, low-pressure cuff, subglottic suction line, and minimally traumatic insertion system



REF 424 

TRACOE experc Set vario XL including REF 461-P extra-long, with adjustable neck flange, low-pressure cuff and minimally traumatic insertion system

To order, please specify: **REF + Size**, e.g. **REF 421-08**